


Underwater Agitator for up
to 18 butyrometers

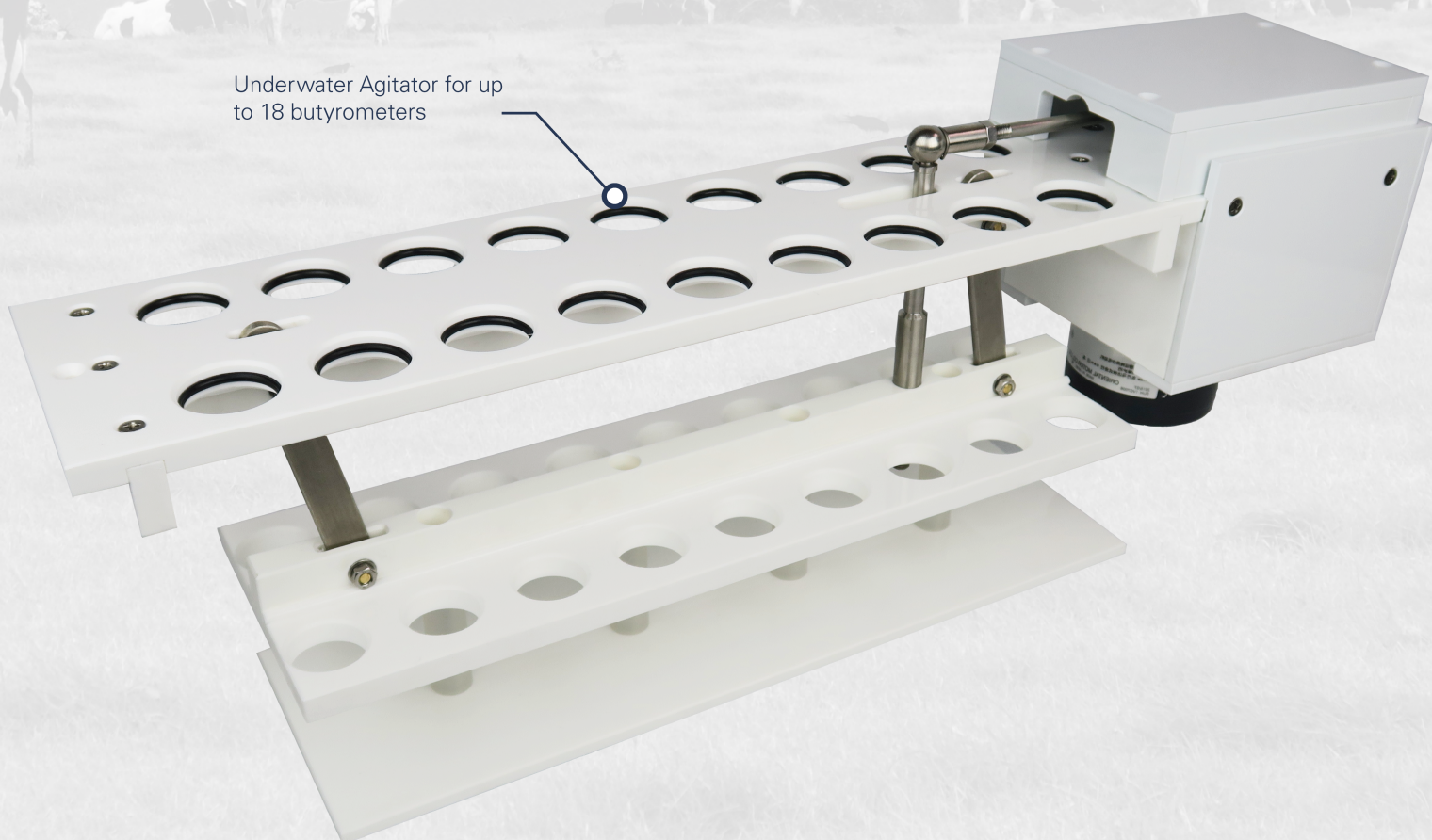

The Aerne Analytic Underwater Agitator is used for an equal intermix of samples in butyrometers at a constant temperature. It works with high precision and safety standards. The agitator can be inserted in waterbaths that have an effective space of 39x25x20cm and 30x15x20cm. The housing

of the Underwater Agitator is made of white polypropylene. It is equipped with an induction motor (shaking 200 UPM). Also in the product range of Aerne Analytic are devices for a very precise handling of butyrometers to analyse fat in laboratories.