


Foil Cutter

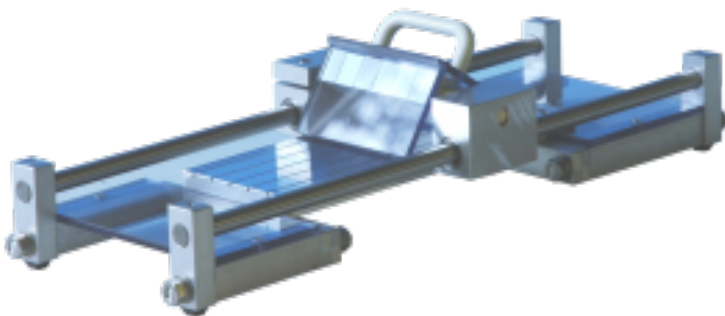

Application:

The Aerne Analytic Foil Cutter is used to cut film strips in exactly the same sizes. The strips are used for quality assurances in the film manufacturings. By means of for example tensile tests the suitability can with the help of a stress-strain diagram be analysed and determined. The cutter works with a ball bearing slide and its hardened cutting blades can be used multiple by reversing them. The Foil Cutter is produced in serial-production to cut four, five or eight films with a width of 15mm and a length of 200mm.

Technical Data:

Cutter with clamping device for foil. Carriage running in ball bearings. Cutting knife can be turned and used four times. Strip length: 200 mm

Order Details:

Part:

Foil Cutter for 4 stripes at 15mm

Foil Cutter for 5 stripes at 15mm

Foil Cutter for 6 stripes at 15mm

Customized Dimensions on Request!

Order-No.

FSG 4

FSG 5

FSG 6

Flyer:



You could also be interested in



Verpackungstester

Aerne Analytic e. K. - Dekan-Schmid-Str. 11- D-89264 Weißenhorn
Tel: +49 (0)7309-92 75 077 -Fax: +49 (0)7309-92 75 076
email: info@aerne-analytic.de

