


Series Waterbath for ex-protection areas

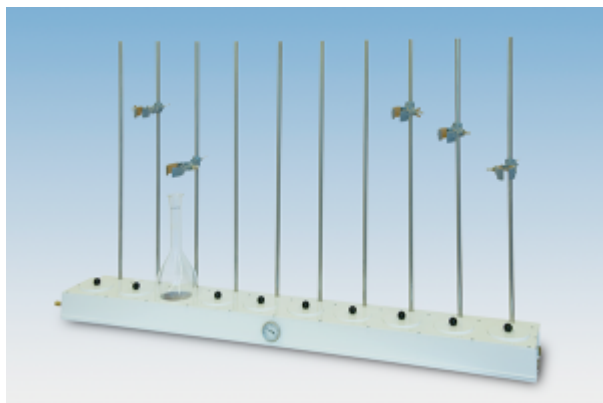

Application:

The Aerne Analytic Series Waterbath for explosion-proof areas are used for extraction or distillation with ether or other inflammable mediums. The heating occurs indirectly by a circulator which is positioned outside of the danger area. The Series Waterbaths of Aerne Analytic are equipped with openings for round or flat-bottomed flasks with removable caps, support rods, clamps and a control dial thermometer. Series Waterbaths can be produced in different sizes and custom specifications.

Technical Data:

Casing: Polypropylene (PP) white Material strength 10 mm

Heating medium: Water

Heating element: nickel-plated copper coils

Scope of supply: Series waterbath with drillings for round or flat bottom flasks, removable Covers, support rods 750 mm long, special round clamps display reference thermometer.

For heating we recommend:

[Revolving Bath Extern Model 1P or Model 2P](#)

Order Details:

Series Water Bath for Round bottom flasks	Length x Width/mm	Order-No.
250 ml 4 Pcs.	505 x 200	RWB 4

6 Pcs.	750 x 200	RWB 6
8 Pcs.	1005 x 200	RWB 8
10 Pcs.	1255 x 200	RWB 10
